



The T-booth – a silent room for hearing testing

- transportable
- permanent rubber seals in all joints – no sticky compound
- built-in silencers in roof and walls
- built-in fan and lamp, both in the roof panel
- large window
- unique cord plug
- frames anodized aluminium
- outside surfaces: vinyl laminated cold rolled steel
- low threshold
- electrically and magnetically shielded

The T-booth has many advantages. Most important is, of course, the excellent sound attenuation, particularly at the critical frequency 500 Hz, which guarantees reliable hearing test results.

Although the T-booth is delivered dismantled, it is simple and quick to assemble.

With the T-booth you are always free to choose if you want to have the window to the left or the right and also on which side you want to hang the door, a good feature the day you want to relocate the T-booth.

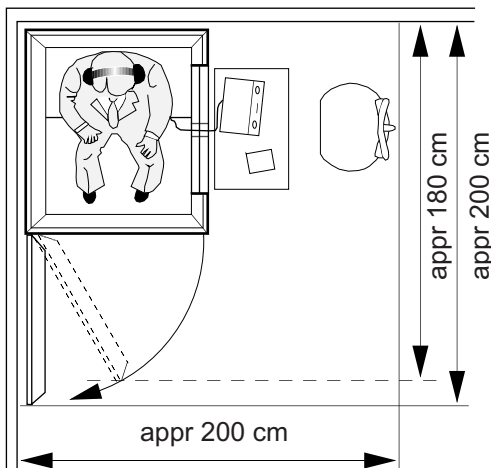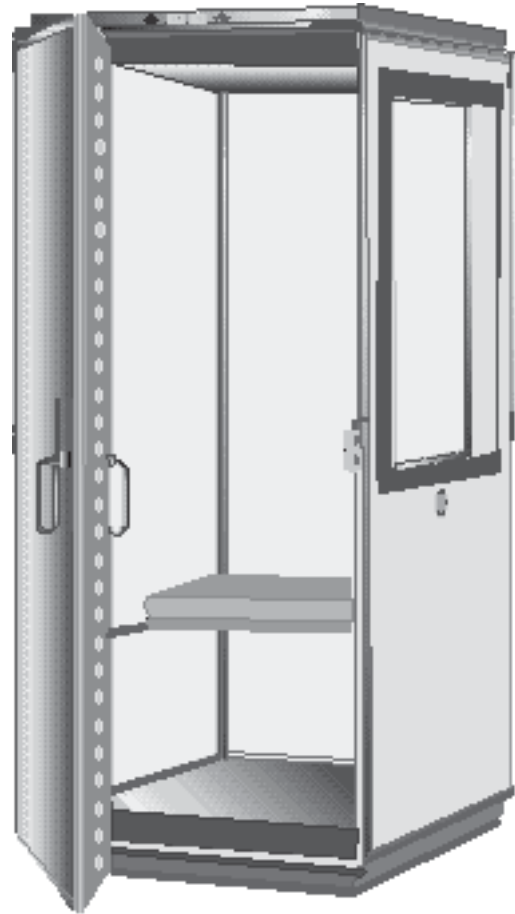

Technical data

DIMENSIONS	WIDTH	DEPTH	HEIGHT
exterior	950 mm	1050 mm	2100 mm
power	220V, 70W		
weight	250 kg net 325 kg (including packing)		

attenuation

Hz	125	250	500	1k	2k	4k	8k
dB	25	33	47	51	51	56	61

test carried out with white noise outside and B&K 2203/613 inside the T-booth